


LB-BME-538KD/LB-BME-538BKGD

Product accessories list

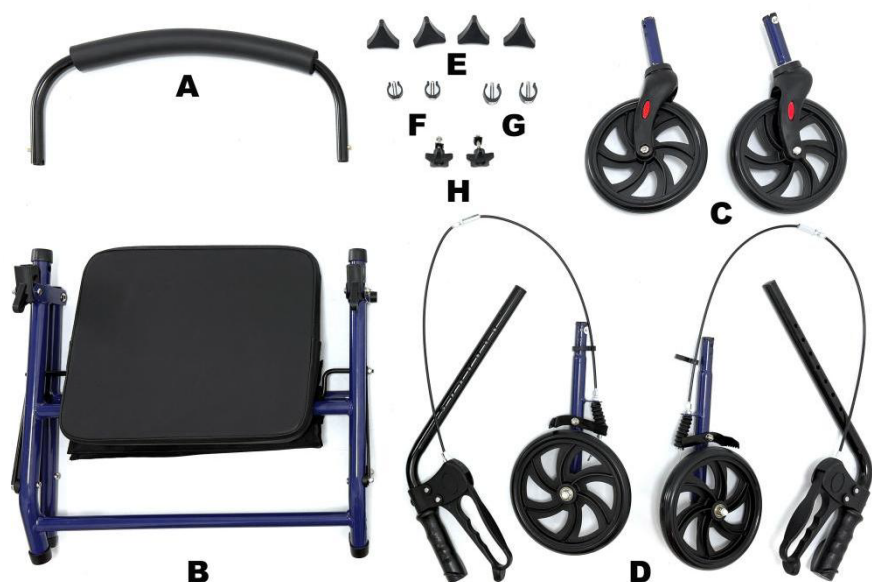

A: Backrest

B: Rollator Frame with Bag

C: Front Wheels*2

D: Handle bars with Brake assemble*2

E: Triangle Screws*4

F: C-pins(S)*2

G: C-pins(L)*2

H: Five Point Kobs with Screw & Curved Washers*2

Important Precautions


Before using your Rollator, you must read and fully understand all of the following safety instructions. This product is a personal mobility aid. Any improper use may cause serious injury to the user or others.

Assembly instructions

Step 1: Open the Rollator



Remove the rollator components from shipping carton. Then, in an upright position, unfold the rollator by pushing down on the folding bracket until the frame is in full extension.

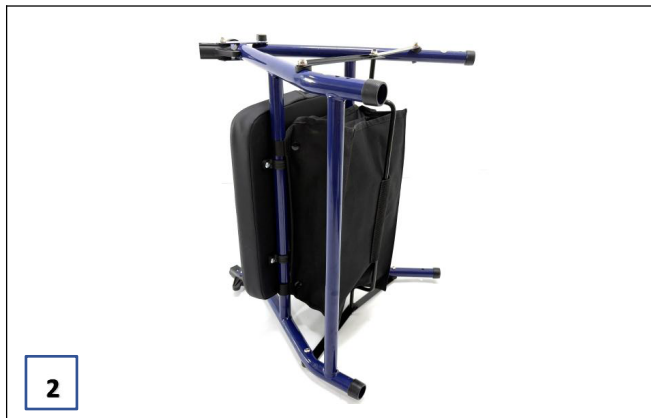

Place the rollator on its side.

NOTE: Be careful when opening the rollator. The moving parts can create pinch points for your fingers.

Step 2: install the wheels

Front wheels



Attach front legs without brakes (C) to the front of the rollator (B). The front legs are marked '1-1', rollator frame '1-1'.



Use C-pins(L) (Part G) to the front legs, then use the triangle screw (Part E) to screw the front legs in place.

NOTE: Note that the C-pins (Part G) should be fully installed and the triangle screws (Part E) should be fully latched.

Rear Wheels & Brakes



With the rollator still on its side, attach rear wheels with brakes (Part D) .

Make sure that the rear wheels and brake shoe are on the outside of the rollator and brake handle (As shown in the image above)



Use the C-pins(S) (Part F) to the rear legs, than use the triangle screw(Part E) to screw the rear legs in place.

NOTE: Note that the c-pins(S) (Part F) should be fully installed and the triangle screws(Part E) should be fully latched.

On the same side insert the brake handle assembly into the frame and set to desired height. Brake cables should be on the outside of the handle bars and frame.



Insert the screw into the hexagonal hole in the frame. Ensure the screw head sits flush within the frame. On the other end of the screw, slide the curved washer on, making sure the curved side faces the frame.



Screw on the four-point knob. Turn the rollator over so the wheels you just installed are on the floor. Then, repeat steps 3-9 on the other side.



Finally install Backrest



Stand the rollator upright so that all four wheels are on the floor. Install the backrest, ensuring the curved portion faces the front of the rollator. Push the backrest into the frame until the locking buttons snap firmly into place. Inspect the backrest to ensure it is secure. Do not use the rollator unless the backrest is securely locked.



Secure the brake cable (As shown in the image above.)



Setup complete. Enjoy!

NOTE: Ensure the handles are set to the same height and that the bag is centered correctly.

Safety Precautions

- 1: Before each use, ensure the rollator is fully opened, the seat is down, and the folding bar under the seat is locked in the down position.
- 2: Check that both height adjustment knobs are fully tightened and set to the same height. Test the rollator's stability before applying your full weight to the hand grips.
- 3: Always test the brakes (loop locks) before use. When the brakes are engaged, the lock bar on the rear wheels should prevent them from turning. If the wheels turn while the brakes are locked, do not use the rollator. Contact your dealer for adjustment. When the brakes are released, the wheels should turn freely.
- 4: Always engage the brakes before sitting. Do not use the seat on an incline or uneven ground.
- 5: Replace tires if they are significantly worn or damaged.

IMPORTANT SAFETY NOTICE

This is a walking aid only. It is NOT to be used as a transportation device. DO NOT allow anyone to push you while you are seated on the rollator.

Folding the Rollator

Find the handle in the bag under the seat.

lift up the handle in the bag.



Brake & Lock Operation


Important: Always apply and release the brakes on both sides at the same time.

1. To Slow Down or Make a Temporary Stop (Service Brake) Gently pull both brake handles upwards towards you. The more you pull, the slower the rollator will move.

2. To Lock the Brakes (Parking Brake)

Press both brake handles down firmly with your palms until you hear and feel them "click" into the locked position. Always check that both brakes are locked before attempting to sit.


3. To Release the Locked Brakes Pull both brake handles upwards until they unlock and return to the neutral, original position.

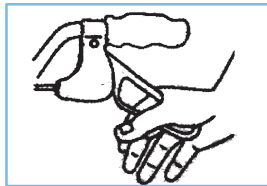

Instructions to adjust the Hand Brake System

Important Safety Information

The brake system installed on this rollator meets all required quality and safety standards. For optimal performance, the brakes have been preset at the factory. However, adjustments may be necessary to suit your personal preference or to compensate for normal wear over time. Regularly check the brake performance and adjust as needed following the instructions below.

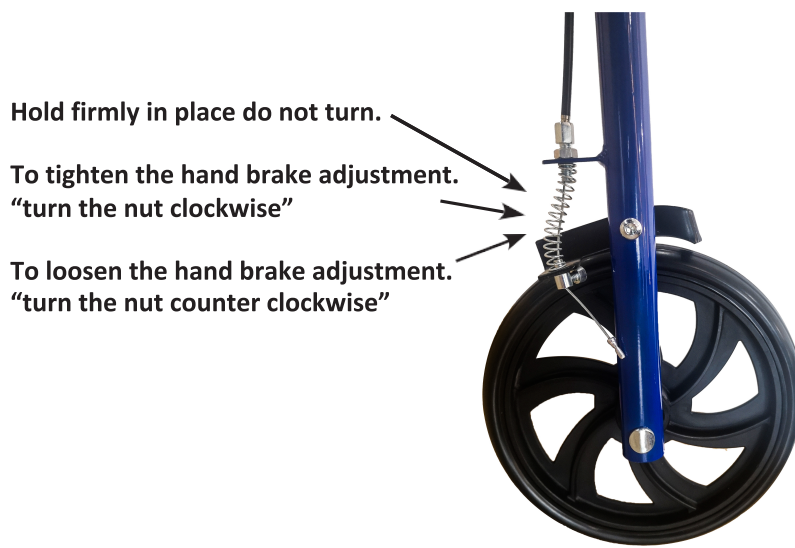

General Safety Instructions

Weight Capacity: DO NOT exceed the maximum user weight capacity or the maximum load capacity for the storage bag as specified.

Not a Transport Device: This product is NOT a wheelchair. NEVER allow anyone to push you while you are seated in the rollator. DO NOT use it to transport people or heavy objects.

Proper Posture: Always walk within the frame of the rollator, not behind it. Keep your body centered between the handgrips to ensure stability.

Both Hands on Grips: Keep both hands on the handgrips at all times when walking to maintain full control and access to the brakes.

Do Not Overload: DO NOT hang heavy items, purses, or bags from the handgrips or the side of the frame. This can cause the rollator to become unstable and tip over.

Be Aware of Surroundings: Watch out for potential hazards on the ground, such as rug edges, electrical cords, thresholds, and spills.

Brake & Seating Safety

ALWAYS LOCK BRAKES BEFORE SITTING: Before you attempt to sit on the seat, you must engage the parking brakes by pushing both brake levers down firmly until they “click” into the locked position. Test that the rollator will not move.

Sit Safely: Turn around and back up to the seat. Use the handgrips for support and slowly lower yourself onto the center of the seat.

Do Not Self-Propel: NEVER attempt to move the rollator by "scooting" with your feet while you are seated.

Maintain Balance While Seated: DO NOT lean forward or to the side to reach for objects while seated, as this can cause a loss of balance and lead to a fall.

Terrain & Environmental Limitations

NEVER Use on Stairs or Escalators: This rollator is strictly forbidden for use on stairs, escalators, or moving walkways.

Use Caution on Slopes: Exercise extreme caution on ramps, inclines, and declines. Use the brakes to control your speed when going downhill.

Avoid Hazardous Surfaces: Be careful on wet, slippery, or icy surfaces. Avoid using the rollator on uneven terrain, loose gravel, thick carpets, or soft ground where wheels can get stuck.

Maintenance & Care

Regular Cleaning: Clean your rollator periodically with a damp cloth and a mild, non-abrasive cleaner. Dry thoroughly.

Check Fasteners: On a monthly basis, check all screws, nuts, and bolts to ensure they remain securely tightened.

Professional Repair: If any part or component—especially brakes and wheels—becomes damaged, worn, or loose, stop using the rollator immediately. Contact your dealer or a qualified professional for repair. Only use manufacturer-approved original replacement parts.

Product Name : Rollator

Model : LB-BME-538KD/LB-BME-538BKGD

Manufacture Date: 202508

Maximum Weight Capacity:300 lbs/136 kg

Brand Owner : Shopper+ INC

Brand Owner Address : 2110 52E AV, LACHINE,
QUEBEC, H8T 2Y3, MONTREAL, CANADA

Email: support@shopperplus.com

Phone: 514-631-5216

Manufacturer :

Manufacturer Address :