LIVING basics®

A brand of Shopper+ Proudly Canadian

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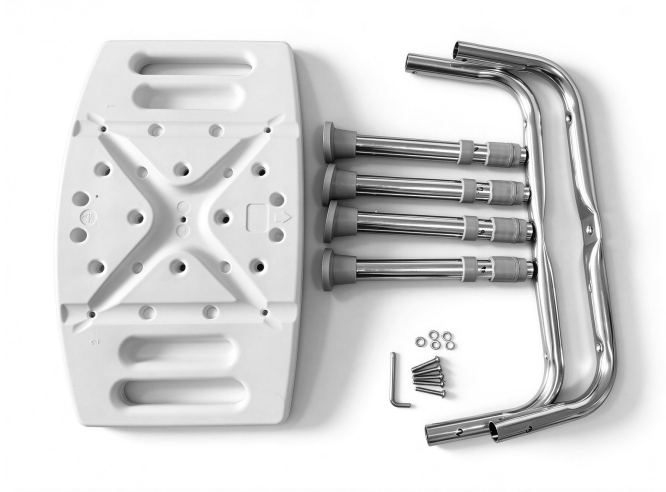
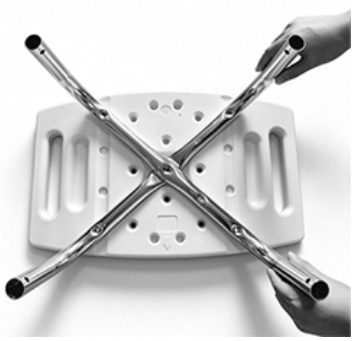
# English

## 1. Important Safety Information

- **Warning:** This chair is for sitting in the bath or shower only. Avoid using it in other environments.
- **Warning:** Ensure the chair is placed on a flat surface to prevent sliding or tipping.
- **Warning:** Avoid excessive force when adjusting the chair's height to prevent damage or injury.
- **Warning:** Regularly check the tightness of all chair parts and adjust if any parts are loose.
- **Warning:** The maximum weight capacity of this chair is 100 kg (220.5 lbs).
- **Warning:** Do not use this chair for elderly or patients who lack balance.

## 2. Installation Instructions

| | |
|---|---|
| **Step 1:**<br>• Remove all parts from the packaging and ensure you have all necessary tools. |  |
| **Step 2:**<br>• Connect the chair legs to the seat base according to the illustrations in the manual. |   |

| | |
|---|---|
| **Step 3:**<br>• Adjust the chair height and ensure all adjustment parts are securely fastened.<br>• Adjust the tips until the chair is parallel to the ground.<br>• Adjust the height of the chair with a push button. |  |
| **Step 4:**<br>• Confirm that all screws and fasteners are tightened, and the chair is stable before use. |  |

# 3. Usage Instructions

- **Preparation:** Ensure the chair is placed securely on a flat surface in the bathroom or shower.
- **Sitting:** Gently sit on the chair and adjust the height to suit your needs.
- **Caution:** Avoid standing or making sudden movements on the chair to prevent slipping.
- **If you are a patient:** Follow your doctor's or therapist's advice when using this chair.

# 4. Intended Use

Intended for use as a shower chair in the bathroom.

- **Indications:**
  - Suitable for users who need to sit while showering, such as elderly individuals or patients who cannot stand for long periods.
- **Contraindications:**
  - Do not use for elderly or patients who lack balance.
  - Not suitable for use outside of the bathroom or shower environment.
  - Not suitable for users exceeding the chair's maximum weight capacity.

- **Notes:**
  - Ensure the floor around the chair is dry to prevent slipping.
  - Regularly check the chair's stability and structural integrity.

## 5. Maintenance and Care

Maintenance and Care

- **Regular Checks:** Periodically check the screws and adjustment parts to ensure there are no loose or damaged components.
- **Cleaning:** Use a mild cleaner and a damp cloth to wipe the seat surface. Avoid corrosive cleaning agents.
- **Storage:** Store the chair in a dry, well-ventilated area after use, avoiding prolonged exposure to moisture.

## 6. Product Specifications

- **Product Name:** Adjustable Height Shower Chair
- **Color:** White
- **Material:** PE Plastic and Corrosion-Resistant Aluminum Frame
- **Model:** LB-BME-2601
- **Maximum Weight Capacity:** 100 kg (220.5 lbs)
- **Seat Adjustment Height:** 17.71 - 21.65 inches (45 - 55 cm)
- **Seat Dimensions:** 16 x 12 inches (40.6 x 30.5 cm)
- **Product Weight:** 1.96 kg (4.32 lbs)

## 7. Warranty Terms

- **Warranty Period:** This product comes with a 90-day warranty. If you find any quality issues or defects during the warranty period, please contact customer service promptly.
- **After-Sales Service:** For warranty service or any inquiries, please contact our customer service team.

## 8. Contact Information

**Brand Owner:** Shopper+ INC

**Brand Owner Address:** 2110 52E AV, LACHINE, QUEBEC, H8T 2Y3, MONTREAL, CANADA

**Email:** support@shopperplus.com

**Phone:** 514-631-5216

**Manufacturer:** FOSHAN DONGFANG MEDICALEQUIPMENT MANUFACTORY(LTD.)

**Manufacturer Address:** NO. 5 XINGFU RD, ZONE A, SHISHAN INDUSTRIAL PARK, NANHAIDISTRICT, FOSHANCITY, GUANGDONG, P.R. CHINA.

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**Fabricant :** FOSHAN DONGFANG MEDICALEQUIPMENT MANUFACTORY(LTD.)